

Mini-Filter Holders

- Suitable for 10 to 20mm filters
- Rod mounted




A mini-filter holder is offered here. It consists of a double clamp mounted on a rod. Filters may be mounted on either side of the rod simply by releasing the spring loaded clamp against the filter. This method is quite secure for filters which rarely need any critical positioning.

Specifications & Tolerances

Dimensions: $\pm 0.2\text{mm}$
 Thickness: 1.5 to 3mm
 Rod length: 20mm
 Material: Aluminum

To hold diameter: 10 to 20mm
 Number of filters: 2
 Finish: Black anodized

Part No. 155-1770

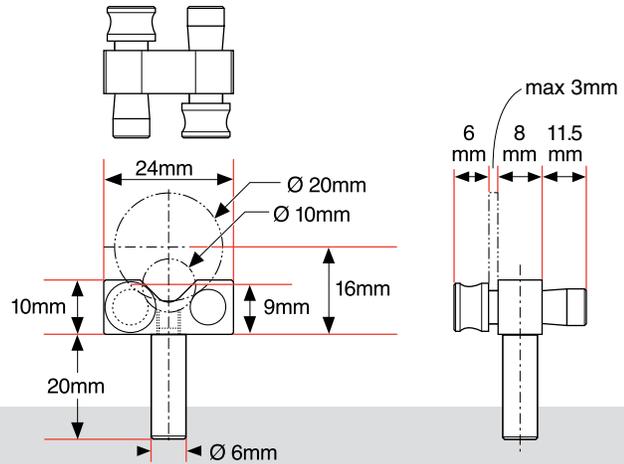

MINI-FILTER HOLDER

ORDERING
&
TECHNICAL SUPPORT
(949) 851-5881
FAX (949) 851-5058
E-MAIL
sales@optosigma.com
WEB
www.optosigma.com

Mini-Filter Holder

To Hold Filter Diameter, D (mm)	Mounting	Price	PART NUMBER
10 - 20	Post		155-1770

